




TARTA® IS STRONG AND FLEXIBLE. BY YOUR SIDE GIVING YOU FREEDOM. IT SAVES YOUR BACK. JUST LIKE A REAL FRIEND.

We were looking for a new idea. We were looking for freedom. We found **Tarta**[®]. An ergonomic back rest, created with the experience of posture professionals with an eye for design and avant-gard materials. Each and every day, we discover the endless possibilities of **Tarta**[®]. For those working in different sectors who share the need to sit for a prolonged period of time, **Tarta**[®] offers the best solution in terms of adaptability. Conceived in a modular mode, with few, simple and sturdy components, **Tarta**[®] back-rests address many needs, from its aesthetic design, adaptability, to the needs of disabled individuals. Everyone can benefit from this innovative solution, now find yours!

It is possible to adjust the inclination of every vertebra with an excursion of 10° , adapting it, for instance, to different degrees of scoliosis. The inclination of every single stave also can be adjusted. The staves being made of aluminum, can be curved at millimetric level for a perfect fit.





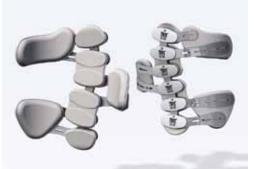




TARTA® HAS INFINITE APPLICATIONS THANKS TO ITS DESIGN AND TECHNICAL MATERIALS. FEEL FREE TO FIND THE BEST ONE FOR YOU

One of the first characteristics of **Tarta**® that will draw your attention is its vast array of solutions and adaptability to suit all your needs. Together with the desired number of vertebrae, other staves of different shape and measure can be applied. On every vertebra, its possible to choose different types of padding. So you can construct your **Tarta**® according to your own individual needs!

All of **Tarta**® materials have been studied and perfected to reach their maximum potential. The harmonic steel used at the core of the vertebrae enables it to retain a 'memory' of the original position if subjected to deformation, giving **Tarta**® flexibility and strength. The more flexible upper part allows for the natural movements of the body, whilst always being able to control the posture. The lower vertebrae, if blocked, offers a valid lumbar support. The padding is fixed to the structure through velcro, so that it can be easily adjusted. The upholstering fabric makes it breathable, giving you total comfort.

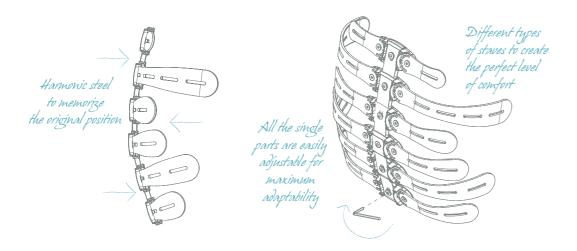

CHAMPIONS CHOOSE THE BEST

There are many sports where a **Tarta**® seat can make the difference. We're testing different solutions to provide our support to the best performance.

EXTRA COMFORT EVERYWHERE

On a couch, when relaxing its possible to gain relief whether you're watching TV or reading a good book.

THE WT AWARD 2011 IN MILAN

We have been awarded a Special Mention in the Accessibility Section from the Province of Milan for our contribution to the improvement of posture for the disabled.














RELAXING IS A QUESTION OF STYLE

There are many places where we sit, often for a very long time. With ${\bf Tarta}^{\otimes}$ you will always be able to maintain the correct posture wherever you sit.

ABOVE AND BEYOND DISABILITY

Easy to mount directly on the back-rest of α wheel chair or pram to find the ideal position.

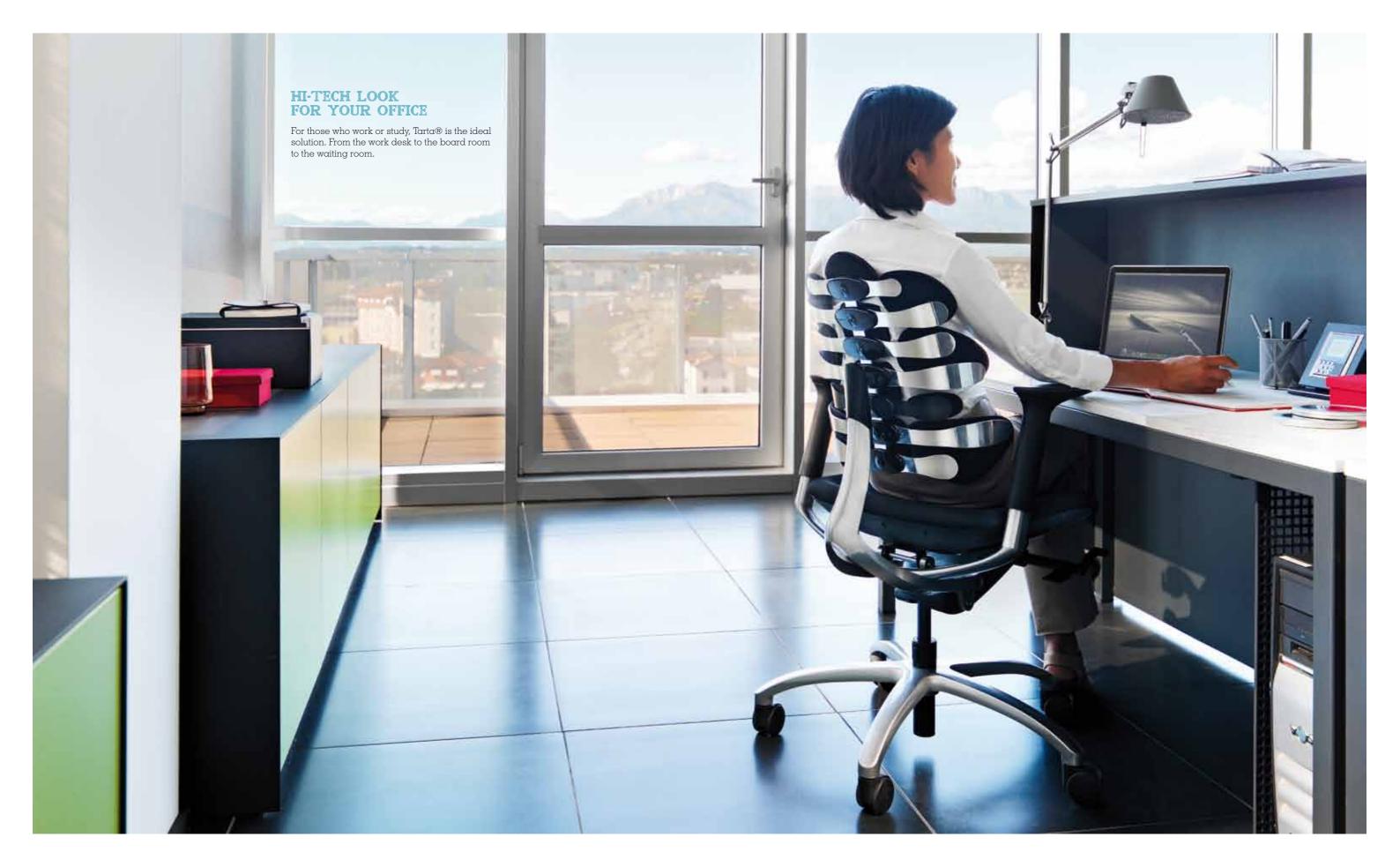

TARTA® EXPRESS THE INNOVATIVE ITALIAN STYLE WITH A SYSTEM OF COMFORTABLE BACK SEATS

The $Tarta^{\circ}$ back-rests reflect a wholly Italian expression of imagination, research and leading design. The advanced types of materials used are given an aesthetically beautiful shape, creating a product that improves and enhances any application. There are many uses for the back-rests, and because we are all different, $Tarta^{\circ}$ offers a wide range of solutions designed to adapt to all of your needs and situations.

Just as with clothing garments, we have created four sizes: S, M, L, XL. In addition to the sizing, there is a choice between having a "Standard" model with semi-embracing slats, giving you a pleasant and comfortable rest; or the "Contour" model inspired by the shaping of racing car seats, offering lateral support thanks to the different lengths and curvature of the slats. Now you can truly find the right solution to fulfill all of your needs.

ARE YOU A COMPANY THAT NEEDS TO CUSTOMIZE THEIR SOLUTION? WE'RE HERE!

We have already operated in cooperation with various companies requiring specific solutions for their items. Some have wanted a certain material, others wanted to change the color, and with others, we have even worked together on creating a particular design to suit their ideas. So, no problem to continue our research with you for any possible customization needs.

STANDARD

CONTOUR

4 VERTEBRAE





5 VERTEBRAE





6 VERTEBRAE







Tarta® Reshapes life

via Cotonificio, 56 33037 Pasian di Prato Udine - Italy T +39 0432 410 089 F +39 0432 722 583

www.tartadesign.it info@tartadesign.it